



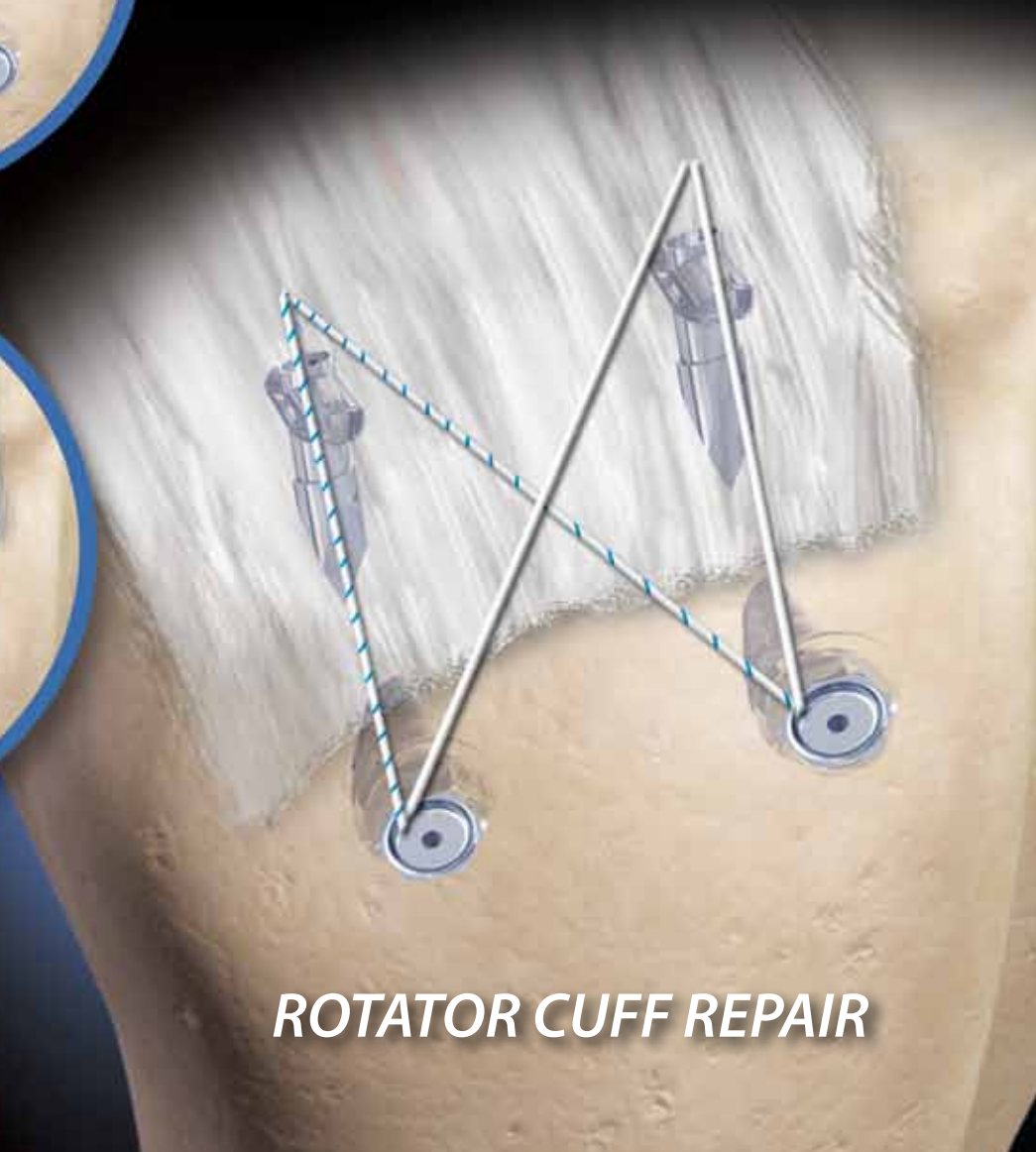

SUTURECROSS[®]

KNOTLESS ANATOMIC FIXATION SYSTEM

- Percutaneous, Transtendon Approach
- Construct Allows for Complete Anatomic Footprint Restoration
- Minimizes Post-Operative Gap Formation
- Arthroscopic or Mini-Open: Double Row, Single Row, PASTA Repairs
- Versatility: All Tears, All Shapes
- Simple, Reproducible Technique

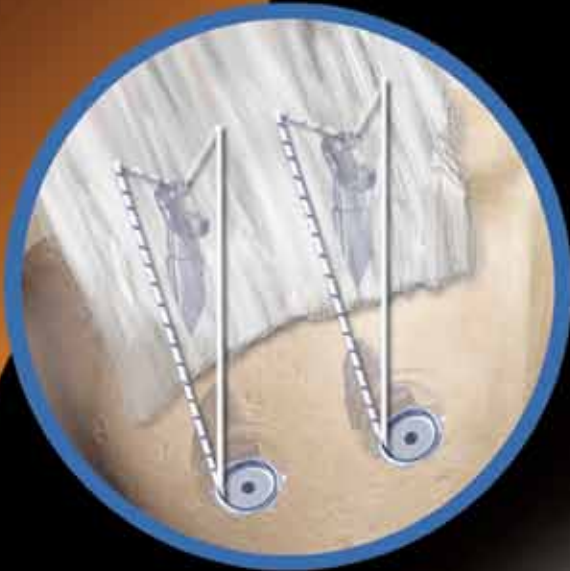

ROTATOR CUFF REPAIR




Product Description:

The KFx® SutureCross® Knotless Fixation System is intended for the fixation of soft tissue to bone using a non-absorbable braided suture. The system consists of a Nail Bone Anchor pre-loaded with a fixed suture in an inserter handle, and a Suture Lock Screw pre-loaded in a delivery handle. All anchors may be placed percutaneously.

Features/Benefits:

- Percutaneous, transtendon or mini-open insertion option.
- Construct allows for complete anatomic footprint restoration.
- Minimizes post-operative gap formation.
- Mini-open or arthroscopic: double row, single row, PASTA repairs.
- Versatility: All Tears, All Shapes.

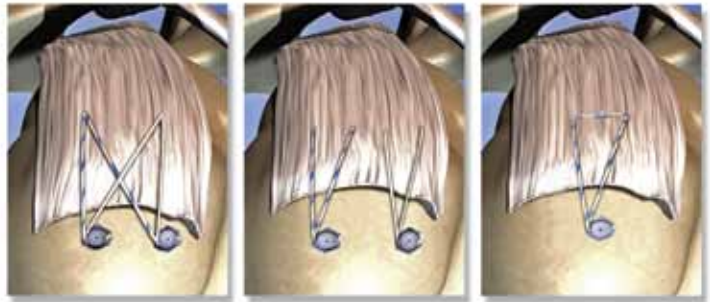

Double Row

Single Row

PASTA



SUTURECROSS® SUTURE NAIL DELIVERY HANDLE



SUTURE LOCK SCREW DELIVERY HANDLE



NAIL BONE ANCHOR



FIXED suture configuration for DOUBLE ROW and SINGLE ROW fixation.

SLIDING suture configuration for PASTA fixations.

SUTURE LOCK SCREW



BONE IMPLANTS:

MODEL NO.	DESCRIPTION
KFX-W-200	One (1), Cobalt Alloy SutureCross® Nail Bone Anchor, One (1), Pre-loaded White #2 Force Fiber® Fixed Suture
KFX-B-200	One (1), Cobalt Alloy SutureCross® Nail Bone Anchor, One (1), Pre-loaded Blue #2 Force Fiber® Fixed Suture
KFX-SW-200	One (1), Cobalt Alloy PASTAFx™ Nail Bone Anchor, One (1), Pre-loaded White #2 Force Fiber® Sliding Suture
KFX-SB-200	One (1), Cobalt Alloy PASTAFx™ Nail Bone Anchor, One (1), Pre-loaded Blue #2 Force Fiber® Sliding Suture

BONE IMPLANTS:

MODEL NO.	DESCRIPTION
KFX-S-155	One (1), Titanium SutureCross® Suture Lock Screw 5.5 mm
KFX-S-163	One (1), Titanium SutureCross® Suture Lock Screw 6.3 mm



5845 Avenida Encinas, Suite 128
 Carlsbad, CA 92008
 Customer Service: 866.883.8718
 760.444.8844 • Fax: 760.444.8855
www.kfxmedical.com

SutureCross® is a registered trademark of KFx Medical, Inc.
 Force Fiber® is a registered trademark of Teleflex, Inc.